

Title: Recording the EKG Lead II

Purpose: The study of the heart's electrical activity is known as electrocardiography. This electrical activity originates from cardiac muscle cells. Electrocardiograms, often known as EKGs or ECGs, are graphs that show how the heart's electrical activity varies.

Procedure:

- 1.To get things started: Before you turn anything on, be sure the IWX/214 unit is plugged in, and that the IWX/214 unit is connected to the laptop by USB cable.Be sure that the CAAMI 504 EEG cable is inserted into the isolated inputs of Channels 1 and 2 of the IWX/214. Be sure that the three color coded lead wires are correctly inserted in the lead pedestal of the CAAMI 504 EEG cable. Insert the connectors on the red, black, and green electrode lead wires into the color coded matching sockets on the lead pedestal of the ECG cable. The white and brown lead wires can be removed and neatly placed in the Iworx case, you do not need them for Lab 10, but after the week is over they will need to be replaced. Once everything is connected, first turn on the laptop and allow it to fully boot up before you turn on the IWX/214 unit. Once the Iworx unit is on, the red indicator light on the Iworx unit should light up and you may hear the USB chime from the laptop if the laptop does not default to mute (many are set to default to mute).
- 2.Open the Labscribe3 program by clicking on the Labscribe3 icon on the desktop. As soon as the program opens, you should see a window pop up that says “Hardware found IWX 214:200824,” click “OK”.
- 3.In the second from the top row (the row that says “File Edit View Tools Settings Advanced External Device Help”), click on the “Settings” tab. About one third of the way down the drop down window should be a tab called “ Human Heart.” Click on that tab and that should lead you to a tab called “ ECG HeartSounds.” Click on that tab and the main window will look like this after you close the pdf file:
- 4.Since Lab 10 is about ECG only, we can hide the lower “Heart Sounds” row by clicking on the symbol to the left of the row label, then clicking on the “ Hide” tab, and then “Yes”.
- 5.Remove the disposable ECG electrodes from its envelope and snap the lead wires onto the electrodes while the electrodes are still on the plastic shield. Instruct the subject to remove all jewelry from their wrists and ankles. Use an alcohol swab to clean a region of skin on the subject’s right wrist, and the inside of both ankles. Let the area dry.

Lab 10

6. Apply the black (1) electrode to the scrubbed area on the right wrist. Repeat Steps 5 and 6 for the inside of the left ankle and the inside of the right ankle, so that the following Lead II is arranged: the black 1 lead is attached to the right wrist, the red +1 lead is connected to the left ankle, the green C or ground lead is connected to the right ankle.

7. Instruct the subject to sit quietly with their hands in their lap. If the subject moves, the ECG trace will move off the top or bottom of the screen. If the subject moves any muscles in the arms or upper body, electromyograms (EMGs) from the muscles will appear on the ECG recording as noise.

8. Click on the Record button, located on the upper right side of the LabScribe Main window. The signal should begin scrolling across the screen. If the ECG appears upside down in Lead II (upside down P, R and T waves), click on the upside down triangle on the far left of "A1:ECG 0.335Hz," then click on the first option "Invert." This should correct the image of your Lead II ECG to be "right side up," but do this ONLY ONCE.

9. When you have a suitable trace, type <Subject's Name> Lead II in the Mark box to the right of the Mark button. Press the Enter key on the keyboard after the recording has started to attach the comment to the data. 10. Click on the AutoScale tab at the upper margin of the ECG channel look for the row that says on the left "A1:ECG 0.335Hz" the AutoScale tab is the second icon after "Hz," it looks like a magnifying glass with a symbol on it. Your recording should look like the figure in step #4. If the ECG waves appear too compressed (too close together), consider clicking the tab above the "Mark" tab that looks like a snow capped pyramid. When the mouse is on top of this tab, it will say "Half Display Time." Clicking this tab will spread out your ECG patterns for step

10. If you overdo that last step, reverse it by clicking on the tab that looks like double pyramids ("Double Display Time") just to the right of the Half Display Time tab.

11. Record for approximately one minute and then click Stop to halt recording. Label one set of the five ECG waves (P, Q, R, S and T). Notice that every cycle is similar but not identical, and the distances between the QRS complexes may alter slightly.




Results:

Lab 10


Time	Lead II	
0	0.00	0.000000
1	0.01	0.000000
2	0.02	0.000000
3	0.03	0.000000
4	0.04	0.000000
...
1454	14.54	0.017423
1455	14.55	-0.026362
1456	14.56	-0.059841
1457	14.57	-0.082315
1458	14.58	-0.097656

1459 rows \times 2 columns

Lab 10

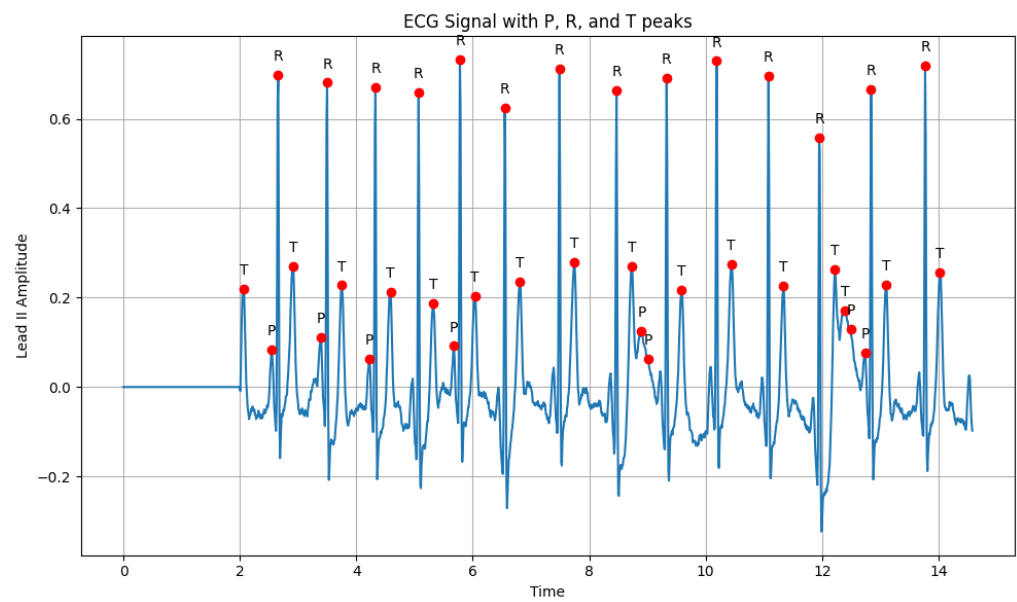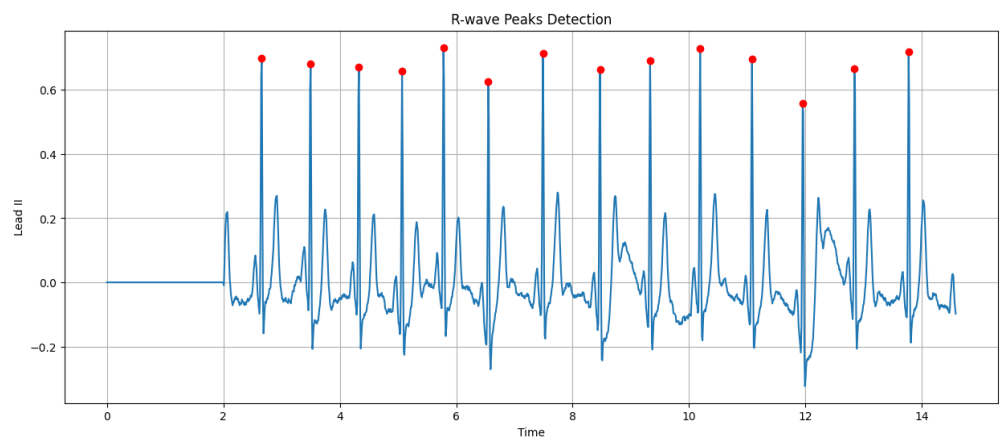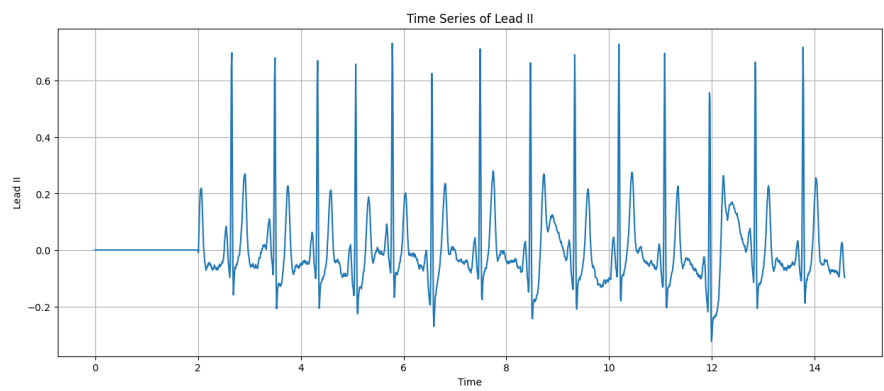

Lab 10

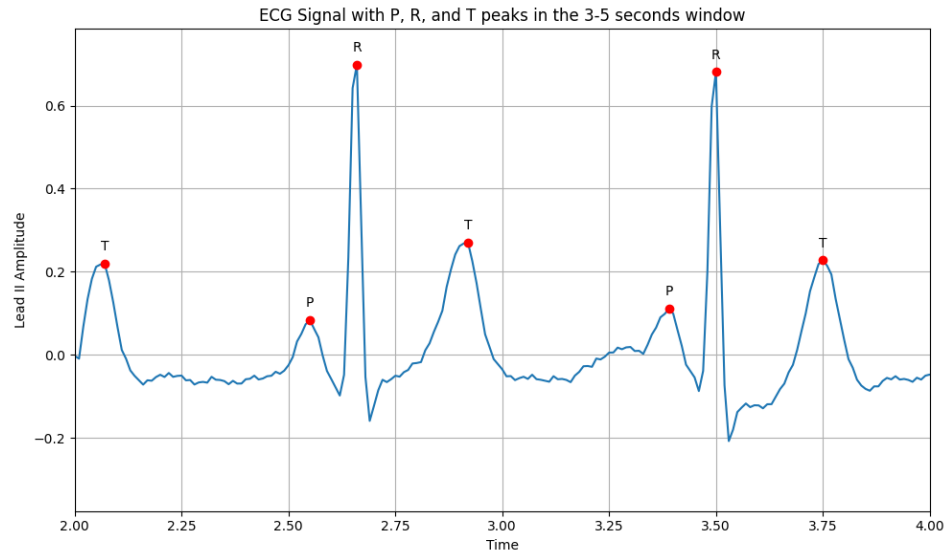

Discussion: The EKG records the electrical signals of the heart. It is a painless procedure to detect abnormalities of the heart.

Conclusion: In conclusion, to be able to identify and explain each component of the EKG is a bit of a tedious task. Also, learning to be able to identify all the instruments used to measure the electrical signals of the heart and to know when there is an abnormality is a great skill to have. It is amazing to see heart rhythms and understand every part of EKG's.